

rIQ 2.0-the Smarter Way to Determine the Refractive Index of Glass Evidence

rIQ 2.0 is the smart way to do glass refractive index measurements that are fast and accurate. Contact CRAIC to learn more.

SAN DIMAS, CA, US, February 11, 2022 /EINPresswire.com/ -- CRAIC Technologies, Inc., the world's leading innovator of UV-visible-NIR microscopy solutions introduce [rIQ 2.0™](#): the smart solution for the analysis of glass trace evidence. rIQ 2.0™, which stands for Refractive Index Quantification, is the result of years of experience in the analysis of trace evidence. rIQ 2.0™ combines sophisticated image analysis software, advanced optical design and electronics to enable criminalists in modern forensic laboratories to measure the refractive index of multiple glass fragments simultaneously, quickly and with the highest accuracy.

"CRAIC Technologies has worked with our customers extensively to build an even smarter way to measure the refractive index of glass quickly and accurately" states Dr. Paul Martin, President of CRAIC Technologies. "Glass fragments are common at crimes scenes. rIQ 2.0™ enables the forensic scientist to measure and compare the refractive index of the smallest fragments of glass with an incredibly high degree of accuracy. And when combined with CRAIC Technologies UV-visible and Raman microspectrophotometers, the refractive index and spectral characteristics of glass fragments can be determined quickly, accurately with the same instrument."

rIQ 2.0™ is an automated system that uses the thermal immersion method, as defined by the standard ASTM E1967, to measure the refractive index of microscopic glass fragments. The system, which incorporates many years of experience with the analysis of glass, allows the user to analyze the refractive index of multiple glass fragments simultaneously and with sophisticated analytical techniques. Statistical analysis methods can also be applied but the instrument is also designed to be user friendly with a short learning curve.




Glass refractive index measurements with rIQ 2.0

riQ 2.0™ is offered as a standalone package, as an add-on package to CRAIC Technologies microspectrophotometers and as an upgrade package for older units already in the field. A standalone package consists of a phase contrast microscope, a high resolution digital camera, the optical interface, a thermal stage, the controlling electronics and the riQ 2.0™ software. The add-on package can be integrated with many CRAIC Technologies microspectrophotometers models, both past and present, to allow them to measure the color, absorbance microspectra™, fluorescence microspectra™, Raman microspectra™ and the refractive index of the smallest of glass fragments.

For more information about riQ 2.0™ for the analysis of the refractive index of glass, visit www.microspectra.com.

Paul Martin
CRAIC Technologies
+1 310-573-8180
sales@microspectra.com

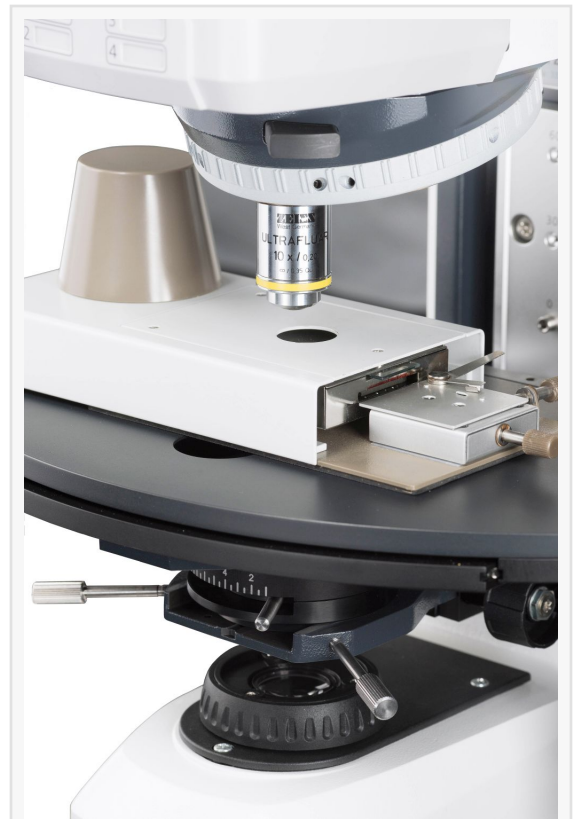

Measuring the refractive index of glass fragments with riQ 2.0

This press release can be viewed online at: <https://www.einpresswire.com/article/562806448>

EIN Presswire's priority is source transparency. We do not allow opaque clients, and our editors try to be careful about weeding out false and misleading content. As a user, if you see something we have missed, please do bring it to our attention. Your help is welcome. EIN Presswire, Everyone's Internet News Presswire™, tries to define some of the boundaries that are reasonable in today's world. Please see our Editorial Guidelines for more information.

© 1995-2022 IPD Group, Inc. All Right Reserved.